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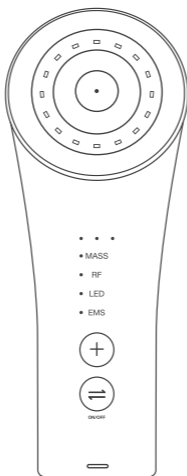
comper

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
Product model: SKC-1S

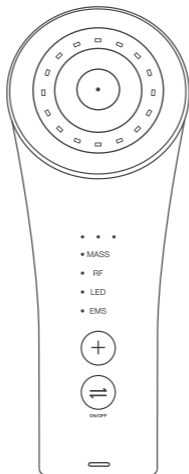comper

To ensure your immediate success, please read this entire User Manual prior to using your device.

## Table of contents

Introduction	2
Product structure, specifications and accessories	4
Before using your device	9
Usage instruction	10
Use Comper Smarkin without App	12
Use Comper Smarkin with App	15
Cleaning, disinfection and storage	20
Frequently asked question	21
Emitted electromagnetic interference	22
Precautions	26
Warranty	27
Contact us	28
Manufacturer	29

## Introduction

### Comper Smarkin Device

Comper Smarkin is a multi-function smart facial machine integrating the RF, EMS, LED and vibration functions that results better anti-aging and face-lifting performance due to their high-efficiency collaboration.

The Smarkin is for anyone wanting to improve skin condition of face and is suitable for most people.

### Intended use

According to the collaboration of functions including RF, LED, EMS and MASS, the skin of face will be adequately stimulated and illuminated except for upper and lower eyelids. Consequently, using Comper Smarkin can improve wrinkles, tighten skin and renew radiance of users'face.

### Clinical application


1. Radio frequency skin tightening (RF) is an aesthetic technique that uses RF energy to heat tissue and stimulate the production of new collagen and elastin in order to remodel tissue and reduce the appearance of fine lines and loose skin.
2. EMS uses low electric currents passing through the skin to produce several therapeutic effects such as muscle toning in the body and micro-lifting of the face.
3. LED is a light source for non-thermal skin rejuvenation because of it can active cells, especially the fibroblasts. LED treatment has been shown to be effective for improving wrinkles and skin laxity.

### Comper Health App

As Smarkin's partner, Comper Health App will analyze your skin health data to help you manage your skin health condition and offers user-defined skin-care routine according to personal skin type and preference. Automatically power off if left idle for 60 seconds.

### Package contents


- Comper Smarkin Device
- Power Adapter
- User Manual

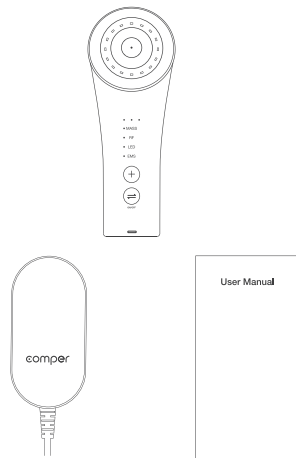
## Product Structure, Specifications and Accessories

### Product structure

- 1.Electrodes
- 2.Temperature Sensor
- 3.LED Therapy
- 4.Power Button / Function-switch Button
- 5.Intensity Up/Down Button
- 6.Intensity Indicator

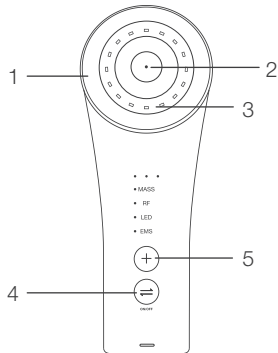- 7.Function Indicator
- 8.Power indicator
- 9.Function Panel: Massage, Radio Frequency, Electro muscle Stimulation, LED Therapy
- 10.DC interface
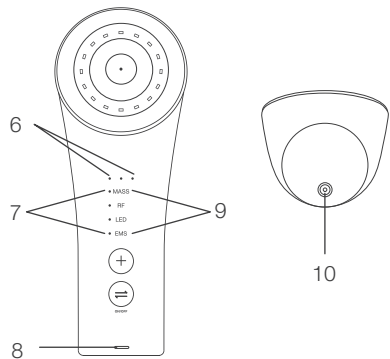
## General specifications

Trade name: Comper Smarkin  
Mode:SKC-1S  
Size: 54\*138\*50(mm)  
Weight: 130±10g  
Input Rated: 100-240V~, 50/60 Hz,  
Output: 12V ~3 A  
Peak Power: 24w  
Functions: MASS, RF, EMS, LED  
Bluetooth:BT4.0,BLE  
Water Resistance:IPX6  
Storage temp:-10° C-55° C  
Working temp:10° C-40° C  
Working Humidity:10%-85%

\* To protect your skin, a temperature sensor is embedded in Comper Smarkin. This sensor ensures the skin's surface will not exceed 42° C / 107.6° F.

## Accessories

HA Essence (not included)


Rub the HA Essence on your treatment area of face before using the Smarkin. This product can offer your skin with sufficient water and help the Smarkin achieve better effect.

Ingredients: Water, butyl glycol, 1, 2 - Pentylene Glycol, poly glutamic acid, beta glucan, methyl gluconic alcohol polyether - 20, hydrolyzed oat protein, sodium hyaluronate, acetyl six peptide - 8, carbomer, hydroxyphenyl propionamide benzoic acid, triethanolamine, phenoxyethanol, ethyl hexyl glycerin, hydroxyl methyl benzene.

## Before using your device

### Hyaluronic Acid Sheet Mask (not included)

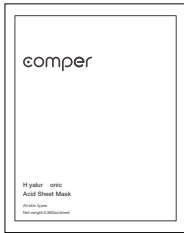Apply the Hyaluronic acid sheet mask after using the Smarkin. The mask can offer your skin with sufficient water and help to improve your skin condition.

Ingredients: Water, glycerin, propylene glycol, butanediol, beta glucan, allantoin, hyaluronic acid, xanthan gum, curacao aloe leaf extract, peony root extract, licorice extract, the hydroxy acetophenone, 1, 2 - diol.

\* Hyaluronic acid sheet mask needs additional purchase.

1. Check the contents of the package, read the user manual comprehensively.
2. Clean up and wipe your face before the operation.
3. Be sure the power adaptor is connected.

## Usage instructions

### Contraindications

- Do not use if you have a history of epilepsy or seizures.
- Do not use if you have had recent facial surgery or other surgical procedure — consult your doctor.
- Do not use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor. If you use a body-worn medical device, such as an insulin pump, consult your doctor prior to use.
- Do not use if you have any active cancer or over any suspicious or cancerous lesion — consult your doctor.
- Do not use if you are pregnant.
- Do not use on infants and children under the age of 18, and keep device and attachments out of the reach of children.
- Do not use if you have metal braces on your teeth.
- Do not use if you have open sores or wounds.
- Do not use if you have prone to rosacea.
- Do not use if you have problematic or sensitive skin.

### Treatment areas


The Smarkin is designed for large surface areas of the face, such as the jawline, cheeks and forehead.

### Caution

Do not use for eyes, eyelids, and surrounding tissue (orbicularis oculi muscle)


Do not use for lips (labium superius oris and labium inferius oris)

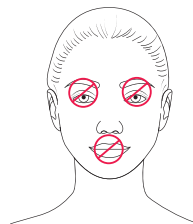
## Use Comper Smarkin without App


1. Rub the HA Essence on the area will be applied by the Comper Smarkin in a moderate quantity and you can add more in following process.

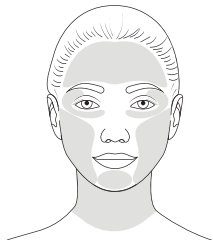2. Being sure the device is connected with the adaptor, hold the power button for 2 seconds to turn on the device, and you can directly start your skincare process.

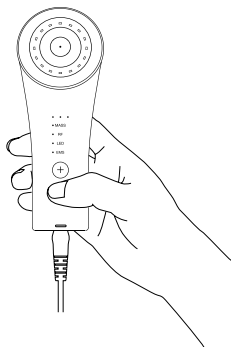
3. Slowly move the device on your treatment areas which are illustrated at "Treatment area" section above.

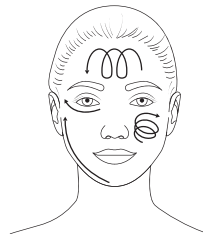
4. The function indicators will indicate the automatic function alteration during the operation, while you can also manually skip the function to the next one by pressing the Function-switch button.

5. Press the intensity up or down button to adjust the level of functional intensity.
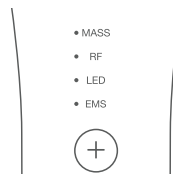
## Use Comper Smarkin with App

6. The Comper Smarkin will stop working as the predefined treatment session being done or time being up and standby for one minute with intensity indicator successively blinking.
7. Hold the Power button for 2 seconds to turn off the device or wait one minute for its automatic turnoff after use.

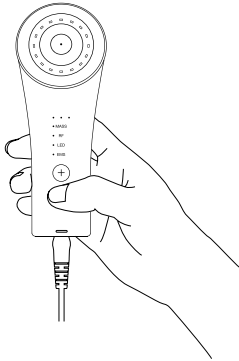
- \* The predefined treatment session usually include three combined functions listed at below consecutively, you will have these normal feelings when you are in each specific part.
- a. MASS+RF function: vibrated and gently heated.
  - b. RF+LED function: gently heated and exposed in red light.
  - c. LED+EMS function: exposed in red light and slightly tingled.

### Bind for first time


1. Search "Comper Health" in App store or scan the QR cord at below, download and open it.

2. Being sure the device is connected with the adaptor and the bluetooth of your smartphone is open.

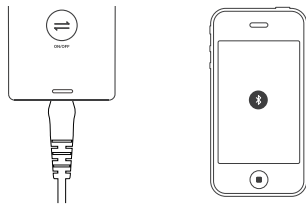
3. Open the App and find "my device", bind the Comper Smarkin device following illustrations below. The App will prompt the success of binding, with the Smarkin buzzing for 1 second and vibrating twice. If binding failed, follow indications on App and rebind it.

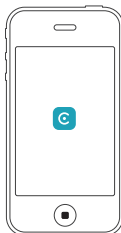4. If you want to control your device with another phone, you need to repeat procedures above again.

### Daily use


1. Be sure the device is connected with the adapter and the bluetooth of your smartphone is open.

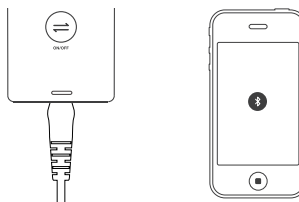2. Open the Comper health App and click "+" at the bottom to start skin care.

3. The Comper Smarkin will stop working and standby for one minute with intensity indicator successively blinking after your skin care routine finished. Hold the Power button for 2 seconds to turn off the device or wait one minute for its automatic turnoff after use.

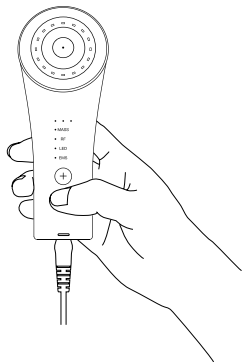
### Restore default skincare mode

The default skincare mode preset in device will be replaced by your personal skincare modes. There are 2 options to quickly go back to default mode.

Option 1 : Directly choose the "Default mode" in Comper Health App.

Option 2: Be sure that the device is connected with the adaptor and the power is turned off. Hold the Power Button and Intensity Up/Down Button at the same time for 3 seconds, the device will vibrate to indicate that the Default Mode is restored.
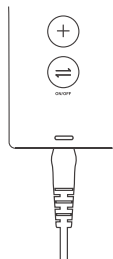
## Cleaning, disinfection and storage.

### Cleaning procedure

1. Turn off the device and disconnect the adapter from it.
2. Wash the device under the water flow with hand or cotton fabric.
3. Wipe the device with a dry cotton fabric or dry it in the air after cleaning.

### Caution

Do not immerse the device into water for long time.

Do not expose the DC interface to water and other liquid.

### Disinfection

Please use cotton cloth with 75% medical alcohol to scrub Smarkin regularly.

### Storage

Do not keep the device in a high temperature or moist place

\*To protect the environment

Do not discard Comper Smarkin with the normal household waste. Dispose of it at an official recycling collection point. Please call your local recycling depot for locations and proper recycling instructions.

## Frequently asked question

### 1.What kind of results will I see, and how quickly will I see them?

Comper Smarkin users have claimed to experience improved facial contour, skin tone and wrinkle reduction. Results and time to onset may vary based on personal skin type and condition. Best results are achieved with regular use. We recommend you to take photos before and after to observe your progress.

### 2.When and how often should I use Comper Smarkin?

There is no limit. We recommend the time before your make-up in the morning and after washing your face in the night. Use the device according to the frequency you can accept and do not overuse.

### 3.Can I use while wearing makeup?

NO. Your skin should be clean and free of any makeup.

### 4.Can I replace HA Essence with other skin care product?

Yes, you can choose other analogous hyaluronic acid product if it could offer similar using experience during device's operation, but Comper is not responsible for any possible problems.

### 5.Is it normal to have a tingling sensation during device's operation?

A slight tingling sensation may occur while using Comper Smarkin. Lowering the intensity, moving the device continuously or increasing the quantity of HA Essence may reduce or eliminate this sensation..

### 6.Can I use comper smarkin when my skin is allergic?

No, Comper Smarkin can't solve the problem of your allergy and even worsen skin's condition.

## Emitted electromagnetic interference

The Comper Smarkin is intended for use in the electromagnetic environment specified below. The customer or the user of the Comper Smarkin should confirm that it is used in such an environment..

Emissions test	Compliance	Electromagnetic environment guidance
RF emissions CISPR 11	Group 1, Class A	The Comper Smarkin uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Table 1: Electromagnetic Environment specified

## Electromagnetic immunity

Emissions test	Compliance	Electromagnetic environment guidance
Electrostatic discharge (ESD) IEC 61000-4-2 6kV contact 8kV air	$\pm 2kV, \pm 4kV, \pm 6kV$ contact $\pm 2kV, \pm 4kV, \pm 8kV$ air discharge	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/ burst IEC 61000-4-4 $\pm 2kV$ for power supply lines	$\pm 2kV$ for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5 $\pm 1kV$ Line to line $\pm 2kV$ Line to earth $\pm 2kV$ Line to earth	$\pm 1kV$ Line to line $\pm 2kV$ Line to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT(>95% dip in UT) for 0,5 cycle 40% UT(60% dip in UT) for 5 cycles 70% UT (30% dip in UT) for 25 cycles <5% UT (>95% dip in UT) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Comper Smarkin requires continued operation during power mains interruptions, it is recommended that the Comper Smarkin be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Conducted RF IEC 61000-4-6 3Vrms 150kHz to 80MHz	3Vrms 150kHz to 80MHz	Portable and mobile RF communications equipment should be used no closer to any part of the Comper Smarkin, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance.
Radiated RF IEC 61000-4-3 3 V/m 80 MHz to 2.5GHz	3 V/m 80 MHz to 2.5GHz	$d = \left[ \frac{35}{f_1} \right] \sqrt{P}$ $d = \left[ \frac{35}{f_1} \right] \sqrt{P_{80MHz \text{ to } 800MHz}}$ $d = \left[ \frac{1}{f_1} \right] \sqrt{P_{80MHz \text{ to } 2.5GHz}}$ <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range B. Interference may occur in the vicinity of equipment marked with the following symbol.</p>
<p>Notes</p> <ul style="list-style-type: none"> <li>At 80 MHz and 800 MHz, the higher frequency range applies</li> </ul> <p>These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p> <p>a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Comper Smarkin is used exceeds the applicable RF compliance level above, the Comper Smarkin should confirm its normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the Comper Smarkin.</p> <p>Over the frequency range 150kHz to 80MHz, field strengths should be less than 3 V/m.</p>		

Table 2: Electromagnetic Immunity Environment

## Recommended safe distance

The Comper Smarkin is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Comper Smarkin can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Comper Smarkin as recommended below, according to the maximum output power of the communications equipment.

Rated maximum out- put power of	Separation distance according to frequency of transmitter (m)		
	$d = \left[ \frac{35}{f_1} \right] \sqrt{P}$	$d = \left[ \frac{35}{f_1} \right] \sqrt{P_{80MHz \text{ to } 800MHz}}$	$d = \left[ \frac{1}{f_1} \right] \sqrt{P_{80MHz \text{ to } 2.5GHz}}$
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.74
0.1	1.17	1.17	2.33
10	3.70	3.70	7.37
100	11.70	1.70	23.30

Table 3: Recommended Safe Distance

## Precautions

- If using Comper Smarkin causes any discomfort, discontinue use immediately.
- If you suspect that you have had a negative reaction from using Comper Smarkin, consult your doctor immediately.
- Please ask your doctor before the use of Comper Smarkin if you have suspected or diagnosed as heart disease , epilepsy, cancer, immune system disease or acute infection.
- The individuals with metal implant or pacemaker around the operational area should avoid of using this device.
- The individuals with allergic or sensation problem on the face should use the device cautiously, or consult the doctor before use.
- Do not use the device over the lips, eyeball, swollen or infected skin, the open wound can't be operated as well except of LED which boots the wound to heal faster.
- The safety of technologies about RF, EMS and LED using for pregnant or postpartum woman have not been proven.
- There may be a perception of flashing lights during the use of your Comper Smarkin, resulting from the stimulation of your optic nerve. If you notice this condition constantly without the use of your Comper Smarkin, consult your doctor.
- Please use the device by keeping away from the equipment with electromagnetic field such as microwave, radio, CT and CCTV.

## Warranty

- 1.Comper promises a one year warranty starting from the date of purchase for free checking and solving all the defects in Comper Smarkin.
- 2.The warranty may not include the damage caused by drop,misuse, abuse ,self-repair and disassemble, natural disasters, improper storage and maintenance.
- 3.Comper provides multiple channels for the global users to find the after-sale service and support. That include the online store, social network , e-mail and official website. Check the "contact us " section in this user manual for more details.
- 4.The warranty do not cover the consequential damages caused by the product, and the any costs in getting the product repaired.

## Contact us

### Type and voice with us in real-time for after sale, technology support and health consultation

1. Search "Wechat" in App store , download and open it .
2. Select "Discover" at bottom , select "Scan QR Code" to scan one of two QR Codes at below:

3. Add Comper official account as your friend, talk with our medical and technical staffs , enjoy the daily communication with Comper's global users in Wechat Group.

### E-mail us for timely feedback within 12 hours

Service@comper.com

### Contact us on social media

Facebook: Comper Healthcare

Instagram: Comperhealthcare

Twitter: Comperhealth

### Get more supports from Comper's website

[www.comper.com](http://www.comper.com)

## Manufacturer

Manufacturer: Comper Chuangxiang (Beijing) Technology Co., Ltd.  
Address: Room 102-103 1st Floor, Building 4, No.1 Kangding Street,  
Daxing District Beijing, China