
User Guide

Comper. SMART BODY FAT SCALE

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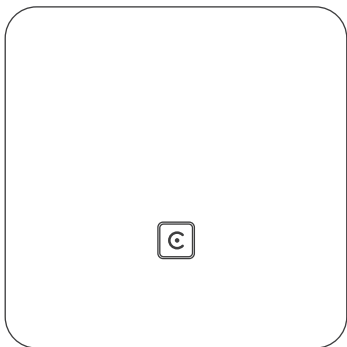Download App

1. Search “Comper Health” in Google Play or App Store, or Scan QR code to download App.


Comper Health


2. Follow the prompts on App to complete registration and binding.
3. To change unit of the weight scale, please enter settings on App.

First Time Use

1. Open the battery cover, put 4 AA batteries to power on the scale

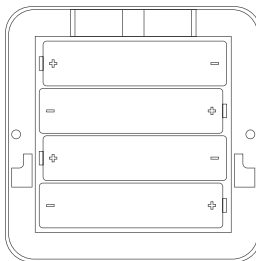
2. Place the scale on the hard, flat surface, step on it to weigh

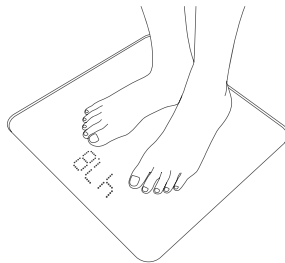
Product Contents


3. 'LO' means you need to change the batteries

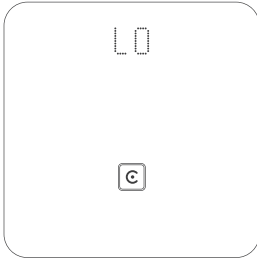
"Err" means the weight exceeds the maximum capability of scale.

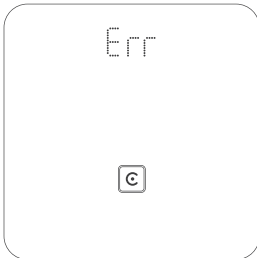One smart body fat scale

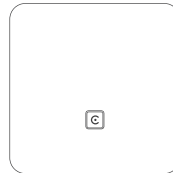Four AAA batteries

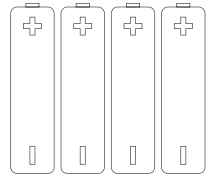One user guide

Measure Range

1. Weigh with bare feet, you can see result of body weight, body fat, body water, muscle mass, BMI (body mass index), BMR (KCAL), bone mass, visceral fat, protein content and body age on Comper Health App.

Note: To ensure accuracy, please be careful not to fit your feet and legs when measuring body fat.

2. Weigh with shoes, you can get body weight and BMI on Comper Health App.

Working Theory

This device is based on the BIA (Bioelectrical impedance Analysis) technology. The reason is that human body fat can not be conductive, but muscle and water are good conductors. The device measures human body bioelectrical impedance to calculate the composition of human body.

Specification

27*80mm Highlight LED Display

Unit: kg / lb

Size: 300*300*30mm

Capacity: 180kg / 400lb

Minimum weight: 5kg / 11 lb

Weight measurement graduation: 0.1kg / 0.2lb

Graduation: Muscle mass,
bone 0.1kg,
Body fat,
body water 0.1%

Age range: 10-99 years old

Height range: 100-220cm

Working humidity: 20%-90%

Storage temperature: -10°C -50°C

Power: DC 6V (4*AAA Batteries)

Caution

Please follow the following warning and reminder to use and maintain products. The company shall not be liable for any damage to persons or articles arising from the operations not in accordance with below contents or the original purpose of the product.

1. Not recommended for individuals with medical implants such as pacemakers;
2. Not recommended for women who are pregnant
3. Make sure your feet are clean before stepping on the scale
4. Do not stand on the edge of the scale or jump on it
5. Do not use the scale in flammable environment, due to possibility of fire or explosion
6. Do not disassemble, destruct, heat or force the scale into a fire, due to possibility of circuit explosion
7. Keep out of children. Adult supervision recommended during child use, to avoid injury
8. Do not put the scale on the wet floor, and ensure feet are dry prior to stepping on the weight panel, due to possibility of slipping
9. Do not place the scale in a wet area, due to possibility of product failure
10. Avoid strong impact, including falls from height, due to possibility of product breakage
11. Do not immerse the scale in water, due to possibility of product failure
12. Please buy the batteries in formal channels

Care and Maintenance

Please use dry and soft cloth to clean the surface of the scale

Wiping the machine by wet cloth with detergent, and clean it with dry cloth

Product Warranty Information

Comper warrants the smart weight scale against defects in quality when used normally in accordance with user guide and ONE (1) YEAR from the date of original retail purchase by the end-user purchase (purchase invoice required):

The following will void the warranty:

1. Product damage caused by negligence of produce specification notice.
2. Self-repair, disassembly, man-made damage, natural disasters and irresistible force.
3. Use more than 1 year.

Legal Disclaimer

Comper Chuangxiang (Beijing) Technology Co. Ltd has the ownership and final interpretation of description of Comper app.

Comper App pursuits for perfect and excellent function, if there is any discrepancy with this specification, subject to the latest version of app.

Comper Chuangxiang (Beijing) Technology Co. Ltd owns the copyright of Comper app.

Copyright

Copyright is owned by Comper (Beijing) Technology Co., Ltd., the Company shall enjoy and reserve all exclusive authority of copyright. No part of the Manual shall be increased or decreased, excerpted, reproduced, altered or imitated without the prior written consent of the Company.

For the content of this Manual, the Company is not responsible for any explicit or implied warranty or guarantee liability. The content change caused by product improvements and technology updates could occur at any time without prior notice. The power of interpretation is vested by the Company.

comper™

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